

ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com Bhubaneswar, Odisha - 751024, INDIA

36-1839: Monoclonal Antibody to CD1a / HTA1 (Mature Langerhans Cells Marker)(O10 + C1A/711)

Clonality: Monoclonal
Clone Name: O10 + C1A/711
Application: FACS,IF,IHC

Reactivity: Human
Gene: CD1A
Gene ID: 909
Uniprot ID: P06126
Format: Purified
Alternative Name: CD1A

Isotype: Mouse IgG1, kappa + Mouse IgG1, kappa

Immunogen Information: Human thymus cells (O10); Recombinant human CD1a protein (C1A/711)

Description

At least five CD1 genes (CD1a, b, c, d, and e) are identified. CD1 proteins have been demonstrated to restrict T cell response to non-peptide lipid and glycolipid antigens and play a role in non-classical antigen presentation. CD1a is a non-polymorphic MHC Class 1 related cell surface glycoprotein, expressed in association with Beta-2 microglobulin. Anti-CD1a labels Langerhans cell histiocytosis (Histiocytosis X), extranodal histiocytic sarcoma, a subset of T-lymphoblastic lymphoma/leukemia, and interdigitating dendritic cell sarcoma of the lymph node. ?When combined with antibodies against TTF-1 and CD5, anti-CD1a is useful in distinguishing between pulmonary and thymic neoplasms since CD1a is consistently expressed in thymic lymphocytes in both typical and atypical thymomas, but only focally in 1/6 of thymic carcinomas and not in lymphocytes in pulmonary neoplasms. Anti-CD1a is reported to be a new marker for perivascular epithelial cell tumor (PEComa).?

Product Info

Amount : 100 μg

Purification: Affinity Chromatography

Content: 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly

toxic.

Storage condition:


Storage condition:

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.

Application Note

Flow Cytometry (1-2ug/million cells); Immunofluorescence (1-2ug/ml); Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 min at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes);

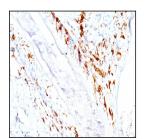

Formalin-fixed, paraffin-embedded human Skin stained with CD1a Monoclonal Antibody (O10+C1A/711).