

ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com

Bhubaneswar, Odisha - 751024, INDIA

36-1964: Monoclonal Antibody to CD63 (Late Endosomes Marker)(NKI/C3 + LAMP3/968)

Clonality: Monoclonal

Clone Name: NKI/C3 + LAMP3/968

Application: FACS,IF,IHC

Reactivity: Human
Gene: CD63
Gene ID: 967
Uniprot ID: P08962
Format: Purified

Alternative Name: CD63,MLA1,TSPAN30

Isotype: Mouse IgG1, kappa + Mouse IgG1, kappa

Immunogen Information: Smooth plasma membrane fraction of MeWo cells (NKI/C3); Recombinant human full-length CD63

protein (LAMP3/968)

Description

This MAb recognizes protein of 26kDa-60kDa, which is identified as CD63. The tetraspanins are integral membrane proteins expressed on cell surface and granular membranes of hematopoietic cells and are components of multi-molecular complexes with specific integrins. The tetraspanin CD63 is a lysosomal membrane glycoprotein that translocates to the plasma membrane after platelet activation. CD63 is expressed on activated platelets, monocytes and macrophages, and is weakly expressed on granulocytes, T cell and B cells. It is located on the basophilic granule membranes and on the plasma membranes of lymphocytes and granulocytes. CD63 is a member of the TM4 superfamily of leukocyte glycoproteins that includes CD9, CD37 and CD53, which contain four transmembrane regions. CD63 may play a role in phagocytic and intracellular lysosome-phagosome fusion events. CD63 deficiency is associated with Hermansky-Pudlak syndrome and is strongly expressed during the early stages of melanoma progression.

Product Info

Amount : 100 μg

Purification: Affinity Chromatography

Content: 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly

toxic.

Storage condition:

Storage condition:

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.


Application Note

Flow Cytometry (1-2ug/million cells); Immunofluorescence (1-2ug/ml); Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes);




ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com


Bhubaneswar, Odisha - 751024, INDIA

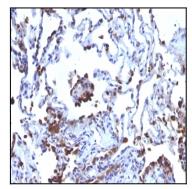

Formalin-fixed, paraffin-embedded human Melanoma metastasized to Lung stained with CD63 Monoclonal Antibody (NKI/C3 + LAMP3/968)

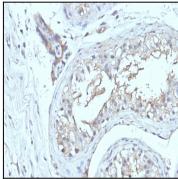

Formalin-fixed, paraffin-embedded human Testicular Carcinoma stained with CD63 Monoclonal Antibody (NKI/C3 + LAMP3/968)

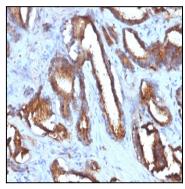

Formalin-fixed, paraffin-embedded human Prostate Carcinoma stained with CD63 Monoclonal Antibody (NKI/C3 + LAMP3/968)