

### 36-11032: Monoclonal Antibody to Clone : IPO-38 (Proliferation Marker)(Clone : IPO-38)

<b>Clonality :</b>	Monoclonal
<b>Clone Name :</b>	IPO-38
<b>Application :</b>	IHC
<b>Reactivity :</b>	Human, Mouse, Rat
<b>Format :</b>	Purified
<b>Isotype :</b>	Mouse IgM, kappa
<b>Immunogen Information :</b>	Spleen cells of a patient with hairy cell leukemia

#### Description


Recognizes a protein of 14-16kDa, which is a novel nuclear antigen of proliferating cells. IPO-38 antigen is present in the nuclei of proliferating cells such as Hodgkin s disease and non-Hodgkin s lymphomas, different forms of leukemias, breast and colorectal carcinomas, and PHA-stimulated lymphocytes. It is not expressed in the cells of non-stimulated lymphocytes and granulocytes. IPO-38 may be a useful marker of cell proliferation during monitoring of tumor progression.

#### Product Info


<b>Amount :</b>	100 µg
<b>Purification :</b>	Affinity Chromatography
<b>Content :</b>	100 µg in 500 µl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly toxic.
<b>Storage condition :</b>	Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid repeated freeze and thaw cycles.

#### Application Note

Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95&degC followed by cooling at RT for 20 minutes);

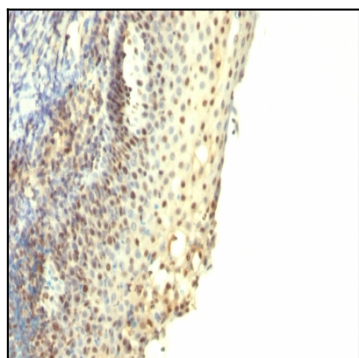

Formalin-fixed, paraffin-embedded human Tonsil stained with IPO-38 Monoclonal Antibody.

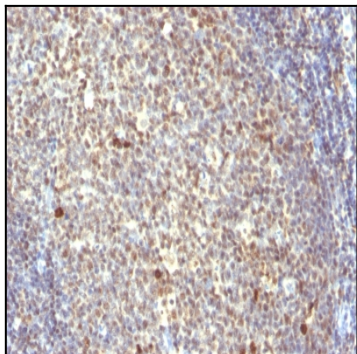

Formalin-fixed, paraffin-embedded human Tonsil stained with IPO-38 Monoclonal Antibody.