

ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com Bhubaneswar, Odisha - 751024, INDIA

36-1506: Monoclonal Antibody to Nucleolin (Marker of Human Cells)(364-5 + NCL/902)

Clonality: Monoclonal

Clone Name :364-5 + NCL/902Application :FACS,IF,WB,IHC

Reactivity: Human
Gene: NCL
Gene ID: 4691
Uniprot ID: P19338
Format: Purified
Alternative Name: NCL

Isotype: Mouse IgG1, kappa + Mouse IgG1, kappa

Immunogen Information: Lysate of SU-DHL-1 Nuclei (364-5); Recombinant human NCL protein (NCL/902)

Description

Recognizes a protein of ~76kDa, which is identified as Nucleolin (NCL). It is the major nucleolar phosphoprotein of growing eukaryotic cells. NCL is located mainly in dense fibrillar regions of the nucleolus. It is found associated with intranucleolar chromatin and pre-ribosomal particles. Human NCL gene consists of 14 exons with 13 introns and spans approximately 11kb. It induces chromatin decondensation by binding to histone H1. It is thought to play a role in pre-rRNA transcription and ribosome assembly.?This MAb can be used to stain the nucleoli in cell or tissue preparations and can be used as a marker of the nucleoli in subcellular fractions. It produces a speckled pattern in the nuclei of cells of normal and malignant cells and may be used to stain the nucleoli of cells in fixed or frozen tissue sections. It can be used with paraformaldehyde fixed frozen tissue or cell preparations and formalin fixed, paraffin-embedded tissue sections.

Product Info

Amount : 100 μg

Purification: Affinity Chromatography

Content: 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly

toxic.

Storage condition:


Storage condition:

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.

Application Note


Flow Cytometry (0.5-1ug/million cells); ,Immunofluorescence (0.5-1ug/ml); Western Blot (0.5-1ug/ml); ,Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT),(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95 °C followed by cooling at RT for 20 minutes),

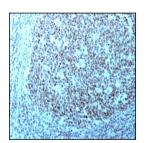


Formalin-fixed, paraffin-embedded human Tonsil stained with Nucleolin Monoclonal Antibody (364-5 + NCL/902




ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com Bhubaneswar, Odisha - 751024, INDIA

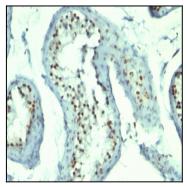


Formalin-fixed, paraffin-embedded human Testicular Carcinoma stained with Nucleolin Monoclonal Antibody (364-5 + NCL/902).

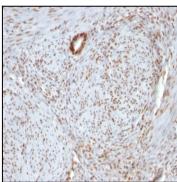

Formalin-fixed, paraffin-embedded human Uterus stained with Nucleolin Monoclonal Antibody (364-5 + NCL/902).

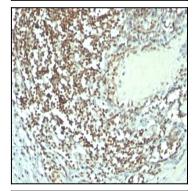

Formalin-fixed, paraffin-embedded human Bladder Carcinoma stained with Nucleolin Monoclonal Antibody (364-5 + NCL/902).

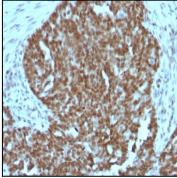


Formalin-fixed, paraffin-embedded human Ovarian Carcinoma stained with Nucleolin Monoclonal Antibody (364-5 + NCL/902).



ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com

Bhubaneswar, Odisha - 751024, INDIA

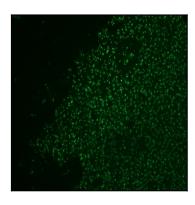

Formalin-fixed, paraffin-embedded human Colon Carcinoma stained with AF488 Conjugate of Nucleolin Monoclonal Antibody (364-5 + NCL/902).