

ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com Bhubaneswar, Odisha - 751024, INDIA

## 36-1657: Monoclonal Antibody to BcI-6 (Follicular Lymphoma Marker)(Clone : BCL6/850)(Discontinued)

Clone Name : Monoclonal
Clone Name : BCL6/850
Application : FACS,IF,IHC

 Gene :
 BCL6

 Gene ID :
 604

 Uniprot ID :
 P41182

 Format :
 Purified

Alternative Name: BCL6,BCL5,LAZ3,ZBTB27,ZNF51

**Isotype:** Mouse IgG1 Kappa

Immunogen Information: Recombinant human bcl-6 protein

## **Description**

Recognizes a protein of 95kDa, which is identified as Bcl-6. Antibody to bcl-6 is helpful in a number of diagnostic settings: (1) In the differential diagnosis of small B-cell lymphoma. Follicular lymphoma will show bcl-6 (and CD10) positivity whereas other small B-cell lymphomas are usually negative. (2) Bcl-6 is an important prognostic marker in diffuse large B-cell lymphomas (DLBCL), where CD10, bcl-6 and MUM1/IRF4 are used to identify germinal center and activated B-cell phenotypes. (3) Bcl-6 can be valuable in distinguishing classical Hodgkin lymphoma from nodular lymphocyte predominant Hodgkin lymphoma (NLPHL). The Reed-Sternberg cells of classical Hodgkin lymphoma are bcl-6 negative whereas the large (?L&H?) cells of NLPHL are bcl-6 positive. In contrast, anti-Bcl-6 rarely stains mantle-cell lymphoma and MALT lymphoma.

## **Product Info**

**Amount** : 100 μg

**Purification:** Affinity Chromatography

Content: 100 µg in 500 µl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly

toxic.

Storage condition:


Storage condition:

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.

## **Application Note**

Flow Cytometry (0.5-1µg/million cells in 0.1ml); Immunofluorescence (1-2µg/ml); Immunohistology (Formalin-fixed) (0.5-1µg/ml for 30 minutes at RT); (Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 10-20 min followed by cooling at RT for 20 minutes); Optimal dilution for a specific application should be determined.

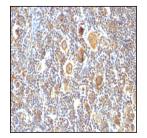


Formalin-fixed, paraffin-embedded human Hodgkin's Lymphoma stained with bcl-6 Monoclonal Antibody (BCL6/850).



ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com

Bhubaneswar, Odisha - 751024, INDIA

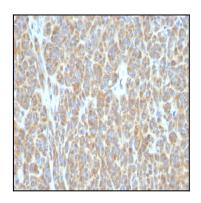

Formalin-fixed, paraffin-embedded human Melanoma stained with bcl-6 Monoclonal Antibody (BCL6/850).