

ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com

Bhubaneswar, Odisha - 751024, INDIA

## 36-1734: Monoclonal Antibody to Transglutaminase II (TGM2)(Clone: TGM2/419)

Clone Name: Monoclonal
Clone Name: TGM2/419
Application: WB,IHC

Reactivity: Human, Mouse, Rat

 Gene :
 TGM2

 Gene ID :
 7052

 Uniprot ID :
 P21980

 Format :
 Purified

 Alternative Name :
 TGM2

**Isotype:** Mouse IgG1, kappa

Immunogen Information: Recombinant full-length human TGM2 protein

## **Description**

Recognizes a 77-85kDa protein, identified as cellular or tissue transglutaminase II (TGase II). Transglutaminases are enzymes that catalyze the crosslinking of proteins by epsilon-gamma glutamyl lysine isopeptide bonds. While the primary structure of transglutaminases is not conserved, they all have the same amino acid sequence at their active sites and their activity is calcium-dependent. The protein encoded by this gene acts as a monomer, is induced by retinoic acid, and appears to be involved in apoptosis. Finally, the encoded protein is the autoantigen implicated in celiac disease. The identification of transglutaminase as the main antigen of endomysium antibodies allows a new diagnostic approach to celiac disease (CD), a genetic, immunologically mediated small bowel enteropathy that causes malabsorption. TGase II is implicated in programmed cell death, signal transduction, drug-resistance, cell growth, endocytosis, insulin secretion, cell adhesion, cataract formation, and wound healing.

## **Product Info**

**Amount :** 100 μg

**Purification:** Affinity Chromatography

Content: 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly

toxic.

Storage condition:

Storage condition:

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.


## **Application Note**

Western Blot (1-2ug/ml); ,Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT),(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95 &degC followed by cooling at RT for 20 minutes),




ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com


Bhubaneswar, Odisha - 751024, INDIA

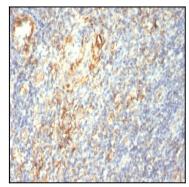


Formalin-fixed, paraffin-embedded human Tonsil stained with Transglutaminase II Monoclonal Antibody (TGM2/419)

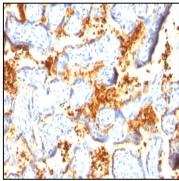

Formalin-fixed, paraffin-embedded human Placenta stained with Transglutaminase II Monoclonal Antibody (TGM2/419)

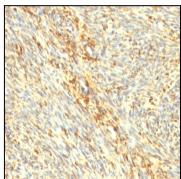

Formalin-fixed, paraffin-embedded human Uterus stained with Transglutaminase II Monoclonal Antibody (TGM2/419)

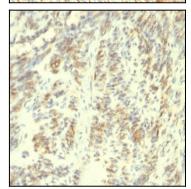

Formalin-fixed, paraffin-embedded human Leiomyosarcoma stained with Transglutaminase II Monoclonal Antibody (TGM2/419)