

ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com

Bhubaneswar, Odisha - 751024, INDIA

## 36-1074: Monoclonal Antibody to Granulocyte-Colony Stimulating Factor (G-CSF)(Clone: SPM468)

Clonality: Monoclonal SPM468 Clone Name: Application: FACS.IF.IHC Reactivity: Human Gene: CSF3 Gene ID: 1440 **Uniprot ID:** P09919 Format: Purified

Alternative Name : CSF3,C17orf33,GCSF Isotype : Mouse IgG1, kappa

Immunogen Information: Nuclei from pokeweed mitogen stimulated human peripheral blood lymphocytes

## **Description**

This MAb recognizes granulocyte-colony stimulating factor (G-CSF) in the cytoplasm of mature granulocytes. It shows no reactivity with any other cell types. Markers of myeloid cells are useful in the identification of different levels of cellular differentiation. It reacts with early precursor and mature forms of myeloid cells. It is useful for the detection of myeloid leukemias and granulocytic sarcomas. It can be used as a marker of granulocytes in normal tissues or inflammatory processes.G-CSF is a pleiotropic cytokine that influences differentiation, proliferation and activation of the neutrophilic granulocyte lineage. The human G-CSF cDNA encodes a 207 amino acid precursor containing a 29 amino acid signal peptide that is proteolytically cleaved to form a 178 amino acid residue mature protein. Two G-CSF s, which are identical except for a three amino acid deletion in the amino-terminus of one form of the protein have been isolated from human cells. Murine and human G-CSF s share 73% sequence identity at the amino acid level.

## **Product Info**

**Amount :** 100 μg

**Purification:** Affinity Chromatography

Content: 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly

toxic.

Storage condition:

Storage condition:

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.


## **Application Note**

Flow Cytometry (1-2ug/million cells); Immunofluorescence (1-2ug/ml); Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT) (Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95&degC followed by cooling at RT for 20 minutes);



ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com

Bhubaneswar, Odisha - 751024, INDIA

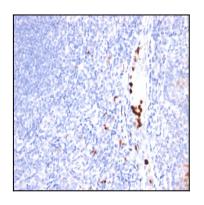

Formalin-fixed, paraffin-embedded human Tonsil stained with G-CSF Monoclonal Antibody (SPM468).