

ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com

Bhubaneswar, Odisha - 751024, INDIA

36-1234: Monoclonal Antibody to HLA-DRB (MHC II)(Clone : SPM288)

Clonality: Monoclonal Clone Name: SPM288

Application: FACS,IF,WB,IHC

 Reactivity :
 Human

 Gene :
 HLA-DRB1

 Gene ID :
 3123

 Uniprot ID :
 P01911

 Format :
 Purified

Alternative Name: HLA-DRB1,HLA-DRB2
Isotype: Mouse IgG2b, kappa

Immunogen Information: Activated human peripheral blood mononuclear cells

Description

This MAb reacts with a 28kDa chain of HLA-DRB1 antigen, a member of MHC class II molecules. It does not cross react with HLA-DP and HLA-DQ. The L243 antibody recognizes a different epitope than the LN3 monoclonal antibody, and these antibodies do not cross-block binding to each other's respective epitopes. HLA-DR is a heterodimeric cell surface glycoprotein comprised of a 36kDa alpha (heavy) chain and a 28kDa beta (light) chain. It is expressed on B-cells, activated T-cells, monocytes/macrophages, dendritic cells and other non-professional APCs. In conjunction with the CD3/TCR complex and CD4 molecules, HLA-DR is critical for efficient peptide presentation to CD4+ T cells. It is an excellent histiocytic marker in paraffin sections producing intense staining. True histiocytic neoplasms are similarly positive. HLA-DR antigens also occur on a variety of epithelial cells and their corresponding neoplastic counterparts. Loss of HLA-DR expression is related to tumor microenvironment and predicts adverse outcome in diffuse large B-cell lymphoma.

Product Info

Amount : 100 μg

Purification : Affinity Chromatography

Content: 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly

toxic.

Storage condition:

Storage condition:

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.


Application Note

Flow Cytometry (1-2ug/million cells); Immunofluorescence (1-2ug/ml); Western Blot (1-2ug/ml); Immunohistochemistry (Formalin-fixed) (0.25-0.5ug/ml for 30 minutes at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes);



ABGENEX Pvt. Ltd., E-5, Infocity, KIIT Post Office, Tel: +91-674-2720712, +91-9437550560 Email: info@abgenex.com

Bhubaneswar, Odisha - 751024, INDIA

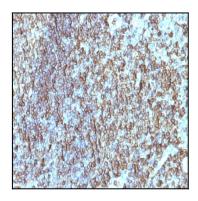

Formalin-fixed, paraffin-embedded human Tonsil stained with HLA-DRB Monoclonal Antibody (SPM288).